

O P E R A T O R Y
F U R N I T U R E



Artizan™ Dental Furniture

The art tradition of furniture design. Continuing the integration of European Design dental furniture into the Midmark product family, Artizan is dental furniture designs crafted of high quality materials that can be customized to meet specific clinic requirements.

Properly designed and manufactured dental casework provides the foundation for efficient operatory function. Artizan dental furniture delivers a distinctive blend of concept and construction specifically designed to meet the unique storage, workflow, infection control and aesthetic requirements of the dental environment.

The Artizan furniture line offers a variety of styles and sizes to suit your needs. If your practice requires a specialized design or size, we welcome the opportunity to work to your specifications. Our experienced representatives are available to assist you with your dental furniture requirements.



TREATMENT STATIONS

Treatment Stations are the primary operatory workstations, providing the needed storage and support at the essential twelve o'clock position. Artizan offers four design series, each with their own set of configuration options.

90 Series

The Artizan 90 Series has been designed for limited space installations. With a full size midsection and upper section combined with a reduced base cabinet depth of less than 17 inches, the 90 Series provides a full function treatment station in a compact package. There is ample storage space for ancillary devices and supplies.

The 90 Series Treatment Station is available in cabinet widths of 36 and 43 inches. The inclusion of a pivoting worksurface offers the assistant increased positioning flexibility and allows for easy left/right conversion.



Monitor Mount
Four monitor mount options are available for the 90 Series. Pivoting model IC-90D shown.



Compact Design
The reduced depth base of the 90 Series allows it to be used in rooms where a full size treatment with pivoting worksurface may not work.



Side Storage
Side storage options include the two-way shelf with tubs. The upper shelf slides forward to present ancillary equipment.

80 Series

Artizan 80 Series Treatment Stations feature a pivoting worksurface, offering increased positioning flexibility for assistant's support. Include a Procenter System as shown for an integrated twelve o'clock instrument delivery system. The dentist's instrument controls are housed behind a bi-fold door in the center base module, enclosed during use but easily accessed when an adjustment is required.

The 80 Series is available in three standard models. The ET3680 is a compact treatment station, 36 inches in width. The ET4380 is a full size treatment station offering ample countertop and storage space. The added depth of the ET4385 provides a side storage compartment for a CPU.

When additional rear cabinetry support is required, the ET4380 offers left or right side expansion modules. The expansion cabinets include an upper storage module, a finished back panel and a sink base module. A drawer module is available as an alternative to the sink base module.



Monitor Mount
The Sliding Track Monitor Mount, model IC-TRX43G, supports the monitor off the counter surface while allowing left/right travel. Additional monitor mount options are available.



Side Kick
The Side Kick side storage option is available with two slide-out shelves or one slide-out shelf and two tubs. The optional Airglide™ opens and closes the side kick with the flip of a switch.



Tub/Tray Support
A tray or tub support and clamp system can be added to the pivoting work surfaces of the 80 and 90 Series for an additional level of instrument support.

70 Series

The Artizan 70 Series Treatment Stations feature a sliding left/right change worksurface. The assistant's worksurface slides out for positioning during procedures. Changing the worksurface left/right position is quick and easy. Lift the small countertop out, slide the worksurface over and reinsert the small top on the opposite side.

The 70 Series is available in three standard models with all the options available in other full size Artizan treatment stations.

The ET3670 delivers the sliding left/right change worksurface in a compact 36 inch width. The ET4370 is a full size treatment station offering ample storage and countertop space. The ET4375 provides dedicated CPU storage, which is the best option when the center base module is required for nitrous rubber goods or other equipment storage.

ET4375
Violin Figured Annigre Cabinet Panels
Canyon Trail Corian® Countertops



Monitor Mounts
The Monitor Arm Mount with Keyboard Tray, model IC-210D, is one of the monitor mount options for the 70 Series Treatment Stations.



Hinged Back Panel
Hinged back panels provide a finished appearance while allowing access.



CPU Side Storage
The ET4375 CPU Side Storage Module provides dedicated CPU storage without sacrificing the front and side storage space.



Easy Left/Right Change in Three Quick Steps...

1. Lift small top upward.



2. Slide tops over.



3. Reinsert top on the opposite side.

60 Series

The Artizan 60 Series Treatment Stations feature a dedicated left or right side sliding worksurface. Although the sliding worksurface does not provide the range of motion offered with a pivoting worksurface, the ease of positioning and smaller room depth requirement can be preferred in many operatory layouts.

The Artizan 60 Series is available in two models; the ET4360 and the ET4365, which includes a CPU side storage module. All the base, midsection, upper and instrument delivery options that are available for other Artizan Treatment Stations are available for the 60 Series.



Side Storage

Five side storage base modules are available for the ET4360, including the Side Kick options which can be upgraded to the Airglide™ system.



Slide-out Worksurface

The sliding worksurface is dedicated to the left or right side and must be specified at time of order.



Slide-out worksurface provides quick and easy positioning.



OPT I O N S



Glass Shelves
Glass shelves maximize the available storage space in the upper module.



Equipment Drawer
A pull-out equipment drawer with a split fold-down front offers storage and convenient access to auxiliary equipment.



CPU Fan
Dedicated CPU storage compartments include a CPU cooling fan.



Two-way Shelves
Side storage two-way shelves pull-out and forward to present support equipment and supplies.



USB Hub
This multiport hub provides easy access for USB connection.




Procenter Control Panel
Handpiece controls are concealed but readily accessible when needed.

CENTRAL STATIONS

Artizan Central Stations are the ideal solution for multi-room plans. Saving valuable floor space, central stations offer storage on both sides for the dentist and assistant. Additional savings are provided as the central station positions one x-ray for utilization in two rooms and one sink for two operatory entries.

Communication between you and your patients is confidential, and HIPAA guidelines require reasonable steps be taken to ensure this confidentiality. The Artizan Tall Boy Central Station addresses these needs by providing the extended height to reduce sound travel and increase privacy in open bay design concepts. And, the optional top trim panel can be modified to accommodate different ceiling heights.

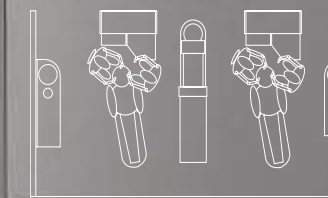

In addition to saving floor space, an open bay layout helps you to economize further. This concept eliminates the need for separate doctor and assistant side sink stations and additional x-ray machines.



Storage Options

Select storage options for each cabinet zone to match your operatory support requirements.



Top Extension

The optional top trim increases the cabinet height, providing an additional degree of privacy.

OPT I O N S



Solid Surface
Corian® solid surface countertop with an undermount stainless steel sink is standard.



Writing Surface
A lift-up writing surface is provided in the sink base storage area.



Waste Chute
Waste chute with easy one hand change bag collar is conveniently located behind the sink.



Dispensing Module
Concealed glove, tissue and mask dispensers are available for the door or sides of the dispenser module.

Zone Options

The Artizan™ Central Station has six designated cabinet zones, each with its own set of options. The sink base comes standard with a Corian® solid surface countertop, an undermount stainless steel sink bowl with matching faucet and soap dispenser. An optional upgrade to an integrated solid surface sink with white faucet and soap dispenser is available. Please refer to our central station ordering sheets for a complete list of options.

SIDE STATIONS

Side station storage modules, trash drops and sinks improve workflow and minimize the need for staff to leave the operatory during a procedure. And, the dentist and assistant can utilize the secondary work surface provided by properly positioned side cabinets.

Artizan Side Stations feature our unique round end storage compartment or round end sink module and are available in seven standard designs. A Tapmaster™ air-less toe kick faucet control is standard, a knee activated control is available. Additional standard features include an undermount stainless steel sink, waste chute and soft close drawers and doors.

Upper storage cabinets and glove/towel dispensers complete the casework suite. The additional cabinetry assists in operatory organization and asepsis by providing a means to keep supplies off countertops and in enclosed storage.



Glove, Towel and Cup Dispenser
Dispenser modules provide a sanitary storage solution for gloves, tissues, masks, towels and cups.



Wall Storage Module
Wall storage modules offer additional operatory storage for bulk items.



CPU Storage Unit
The RS7210 Side Station features a round CPU storage module with cooling fan.



Round End Storage
Functional round end storage modules are available on many models.



O P T I O N S



Round End Storage
Aesthetic, functional round end storage with curved door.



Solid Surface Countertops
Corian® solid surface countertops include an integrated backsplash, a waste chute and undermount stainless steel sink.



Soft Close Drawers
The soft close doors fully extend to reveal the entire contents of the drawers and close smoothly with just a touch.



Waste Chute
A combination waste chute and recessed Sharps container is available.

S T A N D A R D
F E A T U R E S



Concealed, Soft Close Hinges
Concealed, soft close steel hinges provide a clean appearance and allow doors to close smoothly with a slight push.



Tapmaster™
Tapmaster foot activated faucets require no additional utilities and deliver reliable performance.



Stainless Steel Bases
Seamless stainless steel bases are easy to clean and will maintain their good looks through years of use.



Solid Nickel Handles
Non-corroding solid nickel handles provide a distinctive look and are easy to clean and disinfect.



Edge Treatments
Durable color matched 3mm PVC edgebanding and post-formed accent panels.



Drawer Systems
Full extension, soft close drawers extend to reveal the entire contents and close smoothly with just a touch.



Window Etching
Cabinet windows are etched with the Artizan logo. Custom window etching is available.



Solid Surface Countertops
Dupont™ Corian® solid surface countertops are standard. Laminate and DuPont™ Zodiaq® quartz countertops are available.



CPU Storage
Cooling fans are included with CPU storage compartments.

Dental Cabinetry Comparison Chart

★★★★ Excellent ★★★ Fair ★ Poor	Artizan™ Dental Casework	Other Dental Manufacturers	Local Cabinet-Makers
Base Material A good base should be strong enough to provide sufficient support to the cabinet and be made from painted or stainless steel. Seamed and riveted steel can provide strong bases, but may be difficult to clean and an aesthetic compromise. Basic plywood or particle board with laminate covering can delaminate and swell over time. Finally, the base should incorporate adjustable levelers to level cabinets on uneven floors.	★★★ Seamless food grade stainless steel with integrated, adjustable leg levelers	★★★ Powder-coated, painted steel with integrated, adjustable leg levelers -or- ★★ Seamed or riveted steel with integrated, adjustable leg levelers	★ Basic plywood or particle board with laminate covering and wood shims for leveling
Cabinet Frame The frame needs to be strong and durable to withstand the demands of the healthcare environment. Since it will encase many components, it must be built with stability. Modular design allows future integration of new storage options.	★★★ 3/4" thick, 70 - 80 lb. density particle board	★★ 5/8" and 1/2", 45 lb. or greater density industrial grade particle board	★ Low density particle board and various types of plywood
Panel Substrates The substrate is the core material used in a cabinets paneling. Higher density material such as medium density fiberboard (MDF) or industrial grade particleboard is recommended because of its stability and resistance to warping. Other materials such as low density particleboard or plywood are not good options because humidity can cause shrinkage/swelling which will eventually lead to warping, cracking and separation.	★★★ 3/4" thick, 70 - 80 lb. density particle board	★★ 5/8" and 1/2", 45 lb. or greater density industrial grade particle board	★ Low density particle board and various types of plywood
Finishes All exposed surfaces must be covered and sealed to prevent damage from moisture common in the dental environment. Three superior choices are thermally-fused melamine, high pressure laminate and PVC thermofoil (vinyl-clad) covering because they offer the best resistance to wear, staining and moisture. Paints and varnishes are not recommended because they will not stand up to the stringent agents and disinfectants used in a dental office.	★★★ High pressure laminate or thermally-fused melamine	★★★ High pressure laminate or thermally-fused melamine	★★★ High pressure laminate -or- ★ Varnished, painted or left unfinished
Edge Treatments When using high pressure laminates or thermally-fused melamines, an edge treatment will be required to completely seal the panel corners from moisture. Acceptable treatments are 2 or 3 mm thickness PVC, thermofoil vinyl and post formed high pressure laminate. However, all these materials must have adequate equipment for proper application. Watch out for square-edge laminate stripping because it can produce sharp, dangerous corners and when applied with traditional contact adhesive, it can delaminate over time.	★★★ Combination of post form and 3mm PVC	★★ 3 mm and 2 mm PVC, laminate post formed edges and thermofoil coverings	★★ 2 mm PVC (depends on adequate equipment?) -or- ★ Square edge strips
Countertop Options The most popular countertop options are laminates and solid surfaces. Laminates are the most cost-effective while solid surfaces are easy to repair, non-porous and seamless (good for asepsis). Quartz, which is resistant to scratches, chips and cracking and does not require sealing, is making its way into the market.	★★★ Variety of laminates, solid surfaces and quartz options	★★★ Variety of laminates, solid surfaces and quartz options	★★★ Variety of laminates, solid surfaces and quartz options
Type of Assembly The four most popular methods are the glue and dowel, knock down, screw fasteners and staple fasteners. The glue/dowel, the knock down method and mechanical fasteners are the best because they create strong joints. Staple fasteners are unacceptable because joints are weak and the staples are exposed.	★★★ 100% Glue/Dowel assembly	★★★ Glue/Dowel assembly and knock down fasteners	★ Frequently use staple fasteners

Artizan™ Dental Furniture

The art tradition of furniture design. Continuing the integration of European Design dental furniture into the Midmark product family, Artizan is dental furniture designs crafted of high quality materials that can be customized to meet specific clinic requirements.

Properly designed and manufactured dental casework provides the foundation for efficient operatory function. Artizan dental furniture delivers a distinctive blend of concept and construction specifically designed to meet the unique storage, workflow, infection control and aesthetic requirements of the dental environment.

The Artizan furniture line offers a variety of styles and sizes to suit your needs. If your practice requires a specialized design or size, we welcome the opportunity to work to your specifications. Our experienced representatives are available to assist you with your dental furniture requirements.



TREATMENT STATIONS

Treatment Stations are the primary operatory workstations, providing the needed storage and support at the essential twelve o'clock position. Artizan offers four design series, each with their own set of configuration options.

90 Series

The Artizan 90 Series has been designed for limited space installations. With a full size midsection and upper section combined with a reduced base cabinet depth of less than 17 inches, the 90 Series provides a full function treatment station in a compact package. There is ample storage space for ancillary devices and supplies.

The 90 Series Treatment Station is available in cabinet widths of 36 and 43 inches. The inclusion of a pivoting worksurface offers the assistant increased positioning flexibility and allows for easy left/right conversion.



Monitor Mount
Four monitor mount options are available for the 90 Series. Pivoting model IC-90D shown.



Compact Design
The reduced depth base of the 90 Series allows it to be used in rooms where a full size treatment with pivoting worksurface may not work.



Side Storage
Side storage options include the two-way shelf with tubs. The upper shelf slides forward to present ancillary equipment.

80 Series

Artizan 80 Series Treatment Stations feature a pivoting worksurface, offering increased positioning flexibility for assistant's support. Include a Procenter System as shown for an integrated twelve o'clock instrument delivery system. The dentist's instrument controls are housed behind a bi-fold door in the center base module, enclosed during use but easily accessed when an adjustment is required.

The 80 Series is available in three standard models. The ET3680 is a compact treatment station, 36 inches in width. The ET4380 is a full size treatment station offering ample countertop and storage space. The added depth of the ET4385 provides a side storage compartment for a CPU.

When additional rear cabinetry support is required, the ET4380 offers left or right side expansion modules. The expansion cabinets include an upper storage module, a finished back panel and a sink base module. A drawer module is available as an alternative to the sink base module.



Monitor Mount
The Sliding Track Monitor Mount, model IC-TRX43G, supports the monitor off the counter surface while allowing left/right travel. Additional monitor mount options are available.



Side Kick
The Side Kick side storage option is available with two slide-out shelves or one slide-out shelf and two tubs. The optional Airglide™ opens and closes the side kick with the flip of a switch.



Tub/Tray Support
A tray or tub support and clamp system can be added to the pivoting work surfaces of the 80 and 90 Series for an additional level of instrument support.

70 Series

The Artizan 70 Series Treatment Stations feature a sliding left/right change worksurface. The assistant's worksurface slides out for positioning during procedures. Changing the worksurface left/right position is quick and easy. Lift the small countertop out, slide the worksurface over and reinsert the small top on the opposite side.

The 70 Series is available in three standard models with all the options available in other full size Artizan treatment stations.

The ET3670 delivers the sliding left/right change worksurface in a compact 36 inch width. The ET4370 is a full size treatment station offering ample storage and countertop space. The ET4375 provides dedicated CPU storage, which is the best option when the center base module is required for nitrous rubber goods or other equipment storage.

ET4375
Violin Figured Annigre Cabinet Panels
Canyon Trail Corian® Countertops



Monitor Mounts
The Monitor Arm Mount with Keyboard Tray, model IC-210D, is one of the monitor mount options for the 70 Series Treatment Stations.



Hinged Back Panel
Hinged back panels provide a finished appearance while allowing access.



CPU Side Storage
The ET4375 CPU Side Storage Module provides dedicated CPU storage without sacrificing the front and side storage space.



Easy Left/Right Change in Three Quick Steps...

1. Lift small top upward.



2. Slide tops over.



3. Reinsert top on the opposite side.

60 Series

The Artizan 60 Series Treatment Stations feature a dedicated left or right side sliding worksurface. Although the sliding worksurface does not provide the range of motion offered with a pivoting worksurface, the ease of positioning and smaller room depth requirement can be preferred in many operatory layouts.

The Artizan 60 Series is available in two models; the ET4360 and the ET4365, which includes a CPU side storage module. All the base, midsection, upper and instrument delivery options that are available for other Artizan Treatment Stations are available for the 60 Series.



Side Storage
Five side storage base modules are available for the ET4360, including the Side Kick options which can be upgraded to the Airglide™ system.



Slide-out Worksurface

The sliding worksurface is dedicated to the left or right side and must be specified at time of order.



Slide-out worksurface provides quick and easy positioning.



ET4360
Pebble Grey Cabinet Panels
Sage Graystone Corian® Countertops

OPT I O N S



Glass Shelves
Glass shelves maximize the available storage space in the upper module.



Equipment Drawer
A pull-out equipment drawer with a split fold-down front offers storage and convenient access to auxiliary equipment.



CPU Fan
Dedicated CPU storage compartments include a CPU cooling fan.



Two-way Shelves
Side storage two-way shelves pull-out and forward to present support equipment and supplies.



USB Hub
This multiport hub provides easy access for USB connection.



Procenter Control Panel
Handpiece controls are concealed but readily accessible when needed.

CENTRAL STATIONS

Artizan Central Stations are the ideal solution for multi-room plans. Saving valuable floor space, central stations offer storage on both sides for the dentist and assistant. Additional savings are provided as the central station positions one x-ray for utilization in two rooms and one sink for two operatory entries.

Communication between you and your patients is confidential, and HIPAA guidelines require reasonable steps be taken to ensure this confidentiality. The Artizan Tall Boy Central Station addresses these needs by providing the extended height to reduce sound travel and increase privacy in open bay design concepts. And, the optional top trim panel can be modified to accommodate different ceiling heights.




Storage Options

Select storage options for each cabinet zone to match your operatory support requirements.



Top Extension

The optional top trim increases the cabinet height, providing an additional degree of privacy.

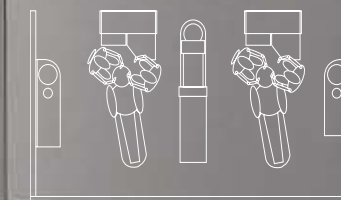

In addition to saving floor space, an open bay layout helps you to economize further. This concept eliminates the need for separate doctor and assistant side sink stations and additional x-ray machines.

OPTIONS



Solid Surface
Corian® solid surface countertop with an undermount stainless steel sink is standard.



Writing Surface
A lift-up writing surface is provided in the sink base storage area.



Waste Chute
Waste chute with easy one hand change bag collar is conveniently located behind the sink.



Dispensing Module
Concealed glove, tissue and mask dispensers are available for the door or sides of the dispenser module.

Zone Options

The Artizan™ Central Station has six designated cabinet zones, each with its own set of options. The sink base comes standard with a Corian® solid surface countertop, an undermount stainless steel sink bowl with matching faucet and soap dispenser. An optional upgrade to an integrated solid surface sink with white faucet and soap dispenser is available. Please refer to our central station ordering sheets for a complete list of options.

SIDE STATIONS

Side station storage modules, trash drops and sinks improve workflow and minimize the need for staff to leave the operatory during a procedure. And, the dentist and assistant can utilize the secondary work surface provided by properly positioned side cabinets.

Artizan Side Stations feature our unique round end storage compartment or round end sink module and are available in seven standard designs. A Tapmaster™ air-less toe kick faucet control is standard, a knee activated control is available. Additional standard features include an undermount stainless steel sink, waste chute and soft close drawers and doors.

Upper storage cabinets and glove/towel dispensers complete the casework suite. The additional cabinetry assists in operatory organization and asepsis by providing a means to keep supplies off countertops and in enclosed storage.



Glove, Towel and Cup Dispenser
Dispenser modules provide a sanitary storage solution for gloves, tissues, masks, towels and cups.



Wall Storage Module
Wall storage modules offer additional operatory storage for bulk items.



CPU Storage Unit
The RS7210 Side Station features a round CPU storage module with cooling fan.



Round End Storage
Functional round end storage modules are available on many models.



O P T I O N S



Round End Storage
Aesthetic, functional round end storage with curved door.



Solid Surface Countertops
Corian® solid surface countertops include an integrated backsplash, a waste chute and undermount stainless steel sink.



Soft Close Drawers
The soft close doors fully extend to reveal the entire contents of the drawers and close smoothly with just a touch.



Waste Chute
A combination waste chute and recessed Sharps container is available.

S T A N D A R D
F E A T U R E S



Concealed, Soft Close Hinges
Concealed, soft close steel hinges provide a clean appearance and allow doors to close smoothly with a slight push.



Tapmaster™
Tapmaster foot activated faucets require no additional utilities and deliver reliable performance.



Stainless Steel Bases
Seamless stainless steel bases are easy to clean and will maintain their good looks through years of use.



Solid Nickel Handles
Non-corroding solid nickel handles provide a distinctive look and are easy to clean and disinfect.



Edge Treatments
Durable color matched 3mm PVC edgebanding and post-formed accent panels.



Drawer Systems
Full extension, soft close drawers extend to reveal the entire contents and close smoothly with just a touch.



Window Etching
Cabinet windows are etched with the Artizan logo. Custom window etching is available.



Solid Surface Countertops
Dupont™ Corian® solid surface countertops are standard. Laminate and DuPont™ Zodiaq® quartz countertops are available.



CPU Storage
Cooling fans are included with CPU storage compartments.

Dental Cabinetry Comparison Chart

★★★★ Excellent ★★★ Fair ★★ Poor	Artizan™ Dental Casework	Other Dental Manufacturers	Local Cabinet-Makers
Base Material A good base should be strong enough to provide sufficient support to the cabinet and be made from painted or stainless steel. Seamed and riveted steel can provide strong bases, but may be difficult to clean and an aesthetic compromise. Basic plywood or particle board with laminate covering can delaminate and swell over time. Finally, the base should incorporate adjustable levelers to level cabinets on uneven floors.	★★★★ Seamless food grade stainless steel with integrated, adjustable leg levelers	★★★★ Powder-coated, painted steel with integrated, adjustable leg levelers -or- ★★ Seamed or riveted steel with integrated, adjustable leg levelers	★ Basic plywood or particle board with laminate covering and wood shims for leveling
Cabinet Frame The frame needs to be strong and durable to withstand the demands of the healthcare environment. Since it will encase many components, it must be built with stability. Modular design allows future integration of new storage options.	★★★★ 3/4" thick, 70 - 80 lb. density particle board	★★ 5/8" and 1/2", 45 lb. or greater density industrial grade particle board	★ Low density particle board and various types of plywood
Panel Substrates The substrate is the core material used in a cabinets paneling. Higher density material such as medium density fiberboard (MDF) or industrial grade particleboard is recommended because of its stability and resistance to warping. Other materials such as low density particleboard or plywood are not good options because humidity can cause shrinkage/swelling which will eventually lead to warping, cracking and separation.	★★★★ 3/4" thick, 70 - 80 lb. density particle board	★★ 5/8" and 1/2", 45 lb. or greater density industrial grade particle board	★ Low density particle board and various types of plywood
Finishes All exposed surfaces must be covered and sealed to prevent damage from moisture common in the dental environment. Three superior choices are thermally-fused melamine, high pressure laminate and PVC thermofoil (vinyl-clad) covering because they offer the best resistance to wear, staining and moisture. Paints and varnishes are not recommended because they will not stand up to the stringent agents and disinfectants used in a dental office.	★★★★ High pressure laminate or thermally-fused melamine	★★★★ High pressure laminate or thermally-fused melamine	★★★★ High pressure laminate -or- ★ Varnished, painted or left unfinished
Edge Treatments When using high pressure laminates or thermally-fused melamines, an edge treatment will be required to completely seal the panel corners from moisture. Acceptable treatments are 2 or 3 mm thickness PVC, thermofoil vinyl and post formed high pressure laminate. However, all these materials must have adequate equipment for proper application. Watch out for square-edge laminate stripping because it can produce sharp, dangerous corners and when applied with traditional contact adhesive, it can delaminate over time.	★★★★ Combination of post form and 3mm PVC	★★ 3 mm and 2 mm PVC, laminate post formed edges and thermofoil coverings	★★ 2 mm PVC (depends on adequate equipment?) -or- ★ Square edge strips
Countertop Options The most popular countertop options are laminates and solid surfaces. Laminates are the most cost-effective while solid surfaces are easy to repair, non-porous and seamless (good for asepsis). Quartz, which is resistant to scratches, chips and cracking and does not require sealing, is making its way into the market.	★★★★ Variety of laminates, solid surfaces and quartz options	★★★★ Variety of laminates, solid surfaces and quartz options	★★★★ Variety of laminates, solid surfaces and quartz options
Type of Assembly The four most popular methods are the glue and dowel, knock down, screw fasteners and staple fasteners. The glue/dowel, the knock down method and mechanical fasteners are the best because they create strong joints. Staple fasteners are unacceptable because joints are weak and the staples are exposed.	★★★★ 100% Glue/Dowel assembly	★★★★ Glue/Dowel assembly and knock down fasteners	★ Frequently use staple fasteners

We have tried our very best to accurately represent product colors on this brochure. However, due to the nature of printed media, the color of actual products may differ slightly from how they appear on your particular printed piece. In order to confirm your color selection, we strongly recommend that you contact Midmark Customer Service at 1-800-MIDMARK to request a sample before placing your order.

CARB 93120.2 Phase 1 Compliant

Midmark is an ISO 13485 and ISO 9001 Certified Company.

For more information or a demonstration, contact your Midmark dealer or Call: 1-800-MIDMARK Fax: 1-877-725-6495
Outside the U.S.A. Call: 1-937-526-3662 Fax: 1-937-526-8214
or visit our website at midmark.com

© 2009 Midmark Corporation
60 Vista Drive P.O. Box 286 Versailles, Ohio 45380-0286
Products subject to improvement changes without notice
Litho in U.S.A. 007-0732-00 Rev. A1 (1/09)



Because we care.